Appareil de massage
facial anti-âge

GUIDE D'UTILISATION FR / EN / ES

FELICITATIONS !

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THANK YOU !

We appreciate your recent purchase of **SPORT-ELEC®AGELESS WONDER™**.

SPORT-ELEC®AGELESS WONDER™ delivers gentle massages to minimize wrinkles and to restore and tone facial muscles for a more youthful appearance.

SPORT-ELEC®AGELESS WONDER™ :

- Improves facial tone
- Rejuvenates the skin
- Evens out skin tone

Please read the manual fully before using this product. If you need further information visit our website www.sport-elec.com.

SAFETY PRECAUTIONS AND WARNINGS

SPORT-ELEC®AGELESS WONDER™ is extremely safe and simple to use. We assume no responsibility for any harm caused by misuse of the product. Please read the entire instruction manual before using the device. **SPORT-ELEC®AGELESS WONDER™** is intended only for cosmetic use on facial tissue for the stimulation of facial muscles as indicated in the instruction manual for personal beauty purposes. All other uses shall be deemed improper.

Use **SPORT-ELEC®AGELESS WONDER™** only on normal, healthy and clean skin.

Do not apply **SPORT-ELEC®AGELESS WONDER™** :

- *Over swollen, infected or inflamed areas of skin or skin eruption.*
- *After experiencing acute trauma or a recent surgery.*
- *Through the head.*
- *Over your chest.*
- *In the case of cancer.*

- *If you are pregnant.*
- *If you are epileptic.*

Never use **SPORT-ELEC**® AGELESS WONDER™:

- *While driving, operating machinery.*
- *While sleeping.*
- *In a humid environment.*
- *In the bath or shower.*

Never immerse the device in water or any liquid.

Keep the device out of the reach of children or people with diminished physical capacities.

Please get your doctor's permission for using **SPORT-ELEC**® AGELESS WONDER™ in conjunction with Botox® or with facial fillers.

Please get your doctor's permission if you have diabetes or if you suffer from muscle or joint problems.


Only use the application wand in accordance with the illustrations in this manual.

Users should start all sessions in a seated position.

Please consult your doctor before using the device if you have heart disease or cardiac infirmities.


The elderly should use the device only under supervision and at low intensity.

LCD DISPLAY

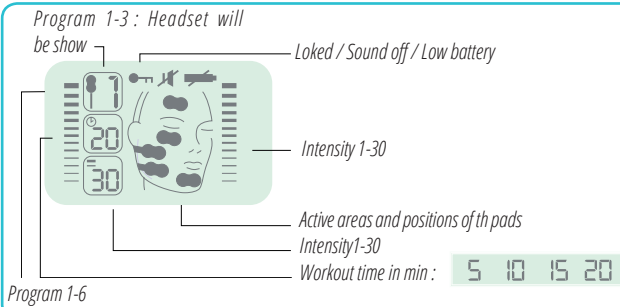HOW TO OPERATE CONTROLS

Indication : Controls :

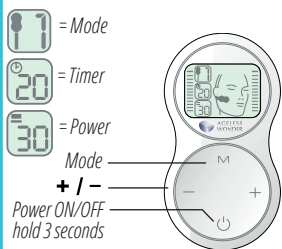
- KEY LOCK TOGGLE ON/OFF :** Press and hold the button together **MODE M** with the “+” while the device wand is operating and you are in intensity mode. The LCD will show the key lock icon.
- MUTE TOGGLE ON/OFF :** aPress and hold the button together **MODE M** etogether with the “-” button for 2 seconds. The LCD will show the mute icon.


Low battery indicator
(SPORT-ELEC® AGELESS WONDER™
requires 2 x AAA batteries)

Before using
SPORT-ELEC® AGELESS WONDER™

1. Always wash your face before use and ensure that it is free from all creams, lotions, oils or make-up.

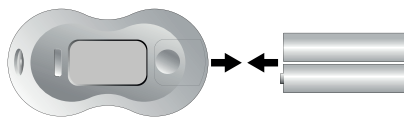2. Open the battery compartment by sliding the cover out. Pull out the tray and place 2 AAA batteries in the proper position. Push the tray back in and replace the battery cover.

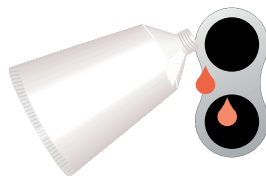3. Use the conductive cream provided and place a small dap of cream on each of the black conductive surfaces. Spread the cream so all the surfaces are covered with a film of cream.


Tip ! For program mode 1 and 2, the correct positioning is directly under your cheek bone towards the back of your face.

SPORT-ELEC® AGELESS WONDER™ is used with a headset or with a handle (optional accessory). We suggest using the headset for program modes 1-3 (follow Option A). For program modes 4-6 use the handle (optional accessory) and refer to option B for installing.

OPTION A/ MODES 1 À 3 Headset use :

Clip the application wands to each end of the headset. Once the application wands are securely installed on the headset, plug the wire into the unit, turn the unit ON and select your options.

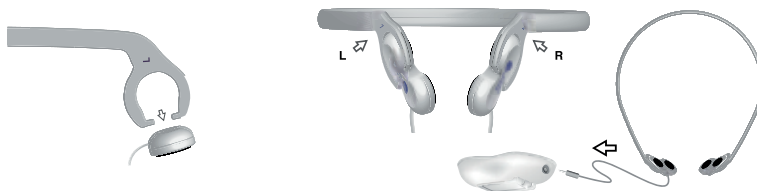OPTION B/ MODES 4 à 5 Handle use :

Clip the optional single application wand into optional handle. Once the application wand is securely installed on the handle, plug the wire into the unit, turn the unit ON and select your options.

PROGRAMS AND APPLICATION AREAS

SPORT-ELEC® AGELESS WONDER™ works on facial muscles that never been stimulated before. Consistency is the key. Use your device for 15 minutes a day for 4 weeks. Use the headset for this stimulation mode. Apply a small dap of cream on each of the black conductive surfaces.

Mode 1 – cheeks– Lift

- *Position the headset horizontally with the application wands below the cheekbone and approx. 1-2 cm in front of the ear.*

This program will act to give you a more lifted appearance and also gives your skin a tighter and smoother look and feel.

Mode 2 – Cheeks – Firm and Tone

- *Position the headset horizontally with the application wands below the cheekbone and approx. 1-2 cm in front of the ear.*

This program will act to give your face a firmer and more toned appearance and also gives your skin a tighter and smoother look and feel.

Mode 3 – Lower Cheek

- *Position the headset horizontally with the application wands below the lower cheekbone and approx. 2-3 cm in front of the ear.*

Select this mode to stimulate the lower cheek, from the corner of the mouth to the upper jaw.

Use the handheld wand for this stimulation mode. Use the wand attached to the cable and plug it into the device. Please note that without the headset, only the application wand with connecting wire will work. We recommend clipping the handle (optional accessory) into the back of the optional single application wand.

Mode 4 – Forehead

- *Use the handheld wand for this stimulation mode.*

Mode 5 - Chin, jaw line

- *Select this mode to stimulate the chin and lower jaw line.*

Mode 6 - Under eye area

- *Select this mode to stimulate the area under your eye.*

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